



This guide is for all Smartphone users.



AbiBird®

How to


Use AbiBird Sensor to Trigger a Request for Contact

Your cared-for person can use the AbiBird sensor to trigger a 'Request for Contact' at any time.

To trigger a 'Request for Contact' alert



Continuously tap or shake the sensor firmly (to activate the vibration sensor) and wait until it beeps and the light turns red.

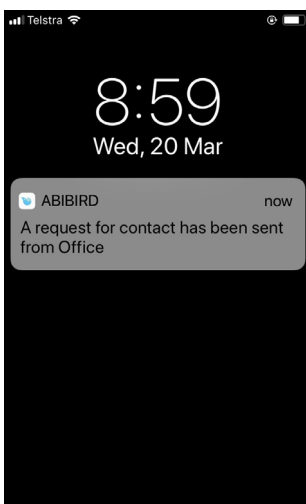


An alert notification and/or email will be sent to the Carer's phone...



This guide is for all Smartphone users.



AbiBird®

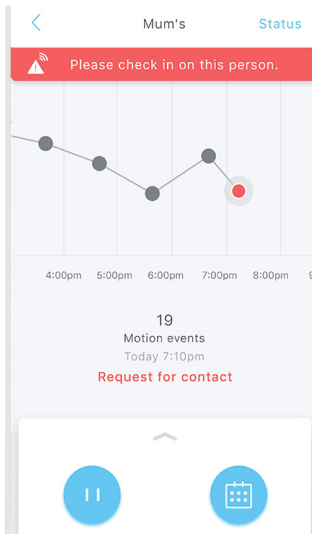

... and the sensor's graph will show a 'Request for Contact' alert (one or more red dots).